


DEQAS July 2018 - 25-OHD Method Means (+/-1SD) for Major Method Groups*

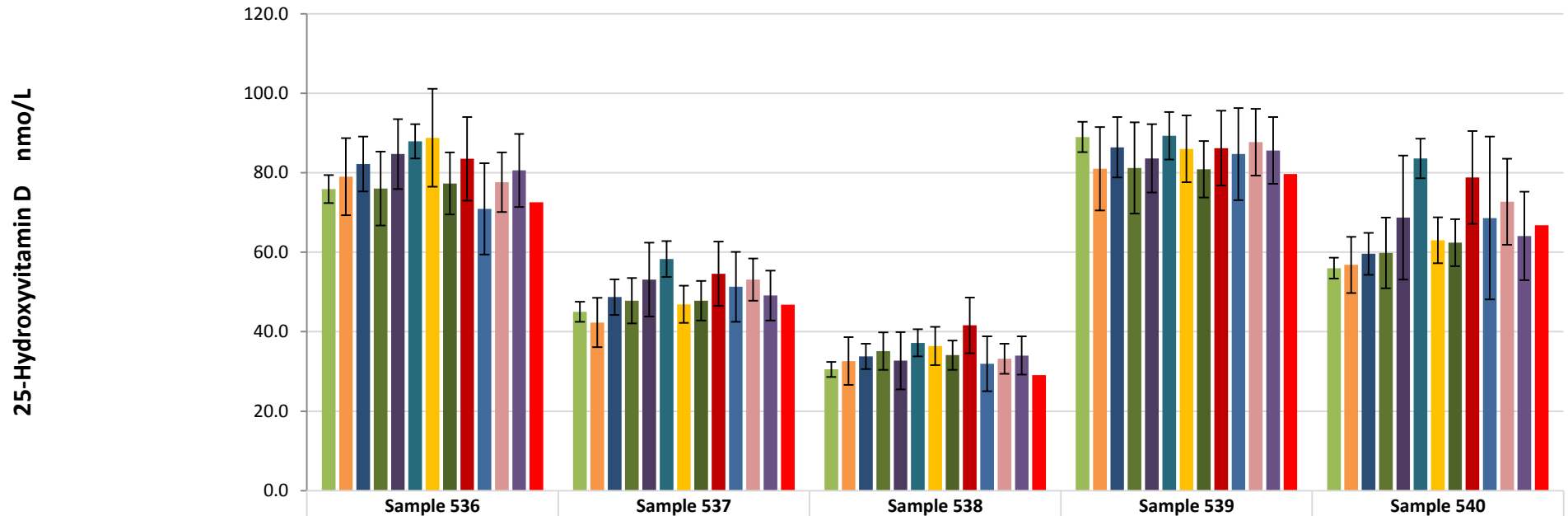

	Sample 536	Sample 537	Sample 538	Sample 539	Sample 540
Abbott Architect New Kit (n=71)	75.9	45.0	30.5	89.0	56.0
Beckman Unicel (n=35)	79.0	42.3	32.6	81.0	56.8
DiaSorin Liaison (n=202)	82.2	48.7	33.8	86.4	59.6
IDS EIA (n=11)	76.0	47.8	35.1	81.2	59.8
IDS iSYS (n=39)	84.7	53.1	32.7	83.6	68.7
IDS-iSYS New (n=19)	87.9	58.3	37.2	89.3	83.6
Roche Total (n=88)	88.8	46.9	36.4	86.0	63.0
Roche Vitamin D total II (n=51)	77.3	47.8	34.1	80.9	62.4
Siemens (n=54)	83.5	54.6	41.6	86.2	78.8
HPLC (n=14)	70.9	51.3	31.9	84.7	68.6
LC-MS/MS (n=141)	77.6	53.1	33.2	87.7	72.7
ALTM (n=775)	80.6	49.1	34.0	85.6	64.1
TARGET VALUE	72.6	46.8	29.1	79.7	66.8

* Methods with 10 or more results returned